

## GLOBAL MEDICAL PRODUCTS

IN CANADA: 1-800-387-6095

AU QUÉBEC: 1-800-655-9445

info@globalmedical.ca



# STAND ON SCALES

## SCALE-TRONIX

5002 Stand On Scale

5202 Stow-a-Weigh Stand On Scale

5125 Portable Stand On Scale

5122 Low Profile Stand On Scale

6102 Flush Mount Floor Scale

# 5002

## Stand-On Scale

- Automatic zero.
- Weight recall.
- Reweigh.
- Cordless.
- 880 lb/400 kg capacity.
- Bright, easy-to-read digital display.

The 5002 is unique. You may turn the 5002 on after the patient is on it and get a correct weight. A paper-tape printer automatically prints out patient weight providing a permanent, accurate record.

The 5002 weighs patients without zeroing in between. The scale automatically returns to zero allowing you to weigh a continuous stream of patients. A "Recall" feature retains the last weight reading after the patient has left the scale—until the next patient is weighed. If the patient gets off the scale and you forget the weight, simply press the "Recall" button. The previous weight is displayed and a "Prior Weight" message blinks. In addition, the press of a button allows reweighing while the patient is still on the scale. An audible signal indicates when a patient's weight is stored in memory.

Our new ergometric design is perfectly balanced to roll effortlessly. The low-profile platform remains motionless while weighing—provides accurate weight on any floor surface, even carpeting. The handrail is ideal for unsteady patients. Patients may hold the handrail without affecting weight readings.

Power cords have been eliminated. Battery life is typically one year, or approximately 20,000 weighings. An optional power supply is available for line cord operation.

### Specifications: 5002 Stand-On Scale

**ACCURACY:** 1/10 LB/100 GM.

**WEIGHING RANGE:** 880 LB/400 KG.

**READOUT:** DIGITAL LED DISPLAY IN POUNDS AND KILOGRAMS.

**POWER SOURCE-CORDLESS:** 6 D-SIZE DISPOSABLE ALKALINE BATTERIES. LINE CORD AVAILABLE.

**AUTOMATIC ZERO:** AUTOMATICALLY RETURNS TO ZERO READY FOR NEXT PATIENT.

**WEIGHT RECALL:** RECALLS LAST WEIGHT.

**REWEIGH:** RECOMPUTES PATIENT'S WEIGHT WHILE PATIENT IS ON SCALE.

**PLATFORM:** LOW PROFILE 18" X 19 1/2" X 1 1/2", WITH HEAVY-DUTY RUBBER WHEELS.

**AUDIBLE SIGNAL:** INDICATES WEIGHT READING IS ENTERED INTO MEMORY.

**LOW BATTERY INDICATOR.**

**WEIGHING MECHANISM:** ELECTRONIC LOAD CELLS.

**SUCCESSIVE WEIGHING:** ALLOWS WEIGHING A CONTINUOUS STREAM OF PATIENTS WITHOUT ZEROING IN BETWEEN.

**WRAP-AROUND HANDRAIL:** HEAVY-DUTY, STAINLESS STEEL.

**HEIGHT GAUGE:** (ACCESSORY). STAINLESS STEEL, IN/CM.

**PRINTER:** (ACCESSORY) BUILT-IN, PAPER TAPE, USES STANDARD THERMAL PAPER.

**COMPUTER CAPABILITY:** RS 232 OUTPUT (OPTIONAL).

### MODEL 5002 - ACCESSORIES:

845010 HEIGHT GAUGE, STAINLESS STEEL, IN/CM.

845010W WALL MOUNT HEIGHT GAUGE. STAINLESS STEEL. IN/CM.

50231 PRINTER

## 5702

### Bariatric Stand On Scale

1000 lbs / 445 Kg weight capacity with a 24" wide x 26" deep platform.

A heavy duty scale with a larger, low profile platform and sturdy handrail.



# 5202

## Stow-A-Weigh® Stand-On Scale

- New, compact fold-up weighing system.
- No bending, no lifting.
- Slim, 4" profile.
- Automatic zero.
- 660 lb/300 kg capacity.
- Cordless.
- Bright, easy-to-read digital display.

This innovative, new Stand-On Scale does a big job while taking up a small space. It's ideal for clinics, triage areas, exam rooms— anywhere that space is limited.

The 5202 mounts snug to the wall providing a slender 4" profile when not in use. Its generous 18" x 18" weighing platform is counter-balanced to fold down easily, no bending or lifting is required.

The digital 5202 automatically turns on and zeros as the platform is opening. Shut-off is also automatic. A 1 1/2" high weighing surface allows easy patient access. Capacity is 660 lb/300 kg. Six D-size alkaline batteries supply approximately 20,000 weighings.

These unique Stow-A-Weigh Scales are accurate to 1/10 lb (100 gm) and, with the touch of a button, will reweigh patients while they're still on the scale. Other advanced features include Weight Recall and Automatic Zero—for weighing a stream of patients without delays. A stainless steel wall-mounted height gauge is optional.

### **Specifications: 5202 Stand-On Scale**

**ACCURACY:** 1/10 LB/100 GM.

**WEIGHING RANGE:** 660 LB/300 KG.

**READOUT:** DIGITAL LED IN LB/KG.

**POWER SOURCE:** CORDLESS, 6 D-SIZE DISPOSABLE ALKALINE BATTERIES.

**AUTOMATIC ZERO:** RETURNS TO ZERO READY FOR NEXT PATIENT.

**WEIGHT RECALL:** RECALLS LAST WEIGHT.

**REWEIGH:** RECOMPUTES WEIGHT WHILE PATIENT IS ON THE SCALE

**PLATFORM:** LOW PROFILE—FOLDS FLUSH AGAINST THE WALL. COUNTER-BALANCED FOR FINGER-TIP OPERATION.

**DIMENSIONS:**

**PLATFORM:** 18"WIDE x 18"DEEP x 1 1/2" HIGH

**OVERALL:** 27" WIDE X 27" HIGH X 4"DEEP.

**LOW BATTERY INDICATOR.**

**WEIGHING MECHANISM:** ELECTRONIC LOAD CELLS.

**SUCCESSIVE WEIGHING:** WEIGH A CONTINUOUS STREAM OF PATIENTS WITHOUT RE-ZEROING.

**HEIGHT GAUGE:** (ACCESSORY) WALL MOUNT, STAINLESS STEEL, IN/CM.

**COMPUTER CAPABILITY:** RS 232 OUTPUT (OPTIONAL).

**MODEL 5202 ACCESSORIES:**

845010W HEIGHT GAUGE (WALL MOUNT), STAINLESS STEEL, IN/CM.



Also available from Global Medical



### **Follett Ice & Water Dispensers**

Healthcare facilities choose Follett ice machines and ice dispensing equipment for long-term reliability, customer-valued features and the industry's best program of service and support.

Now with Agion anti-microbial product protection to minimize biological growth.



### **Metro LifeLine Code Cart**

Designed for a code, not adapted for one. Many of your colleagues have standardized on the LifeLine cart because it allows them the efficiency and control they need in a code situation.



### **Metro Procedure Carts**

Metro family of Procedure carts (Starsys or FlexLine) provide flexible, adaptable mobile storage where you need it. These carts can be configured to suit your need and reconfigured as your needs change. Metro's organized storage systems

# 5125

## Portable Stand-On Scale

- Cordless.
- Low profile with carrying handle.
- Readout contained in platform.
- Lightweight, easy to carry.
- 660 lb/300 kg capacity



The new, ultra-light 5125 Stand-On Scale's compact configuration allows it to be stored easily under a desk or bed. It's equipped with a convenient carrying handle and features a large, low-profile platform for easy access.

The 5125 is designed for user-friendly operation. Just step on the handle, and the scale turns on and zeros.

No bending! The unit automatically returns to zero after each weighing. Patients may be reweighed while still on the scale by just pressing a button. Power for up to 10,000 weighings is supplied by 6 "C" alkaline batteries.

### **Specifications: 5125 Portable Stand-On Scale**

**ACCURACY:** 1/10 LB/100 GM.

**WEIGHING RANGE:** 660 LBS/300 KG. **WEIGHT:** 12 LBS.

**READOUT:** DIGITAL LED DISPLAY IN POUNDS AND KILOGRAMS, SELF CONTAINED IN PLATFORM.

**POWER SOURCE:** CORDLESS, 6 DISPOSABLE ALKALINE "C" BATTERIES FOR 10,000 WEIGHINGS.

**PLATFORM:** 12"w X 17"d X 1 1/4" h.

**WEIGHING MECHANISM:** ELECTRONIC LOAD CELLS.




### **Also Available: 5127 Bariatric Stand-On Scale**

1000 lbs / 445 Kg weight capacity with a 17" x 22" platform

# 5102

## Physician's Stand-On Scale


- Cordless.
- Low profile
- Lightweight, easy to move
- 880 lb/400 kg capacity

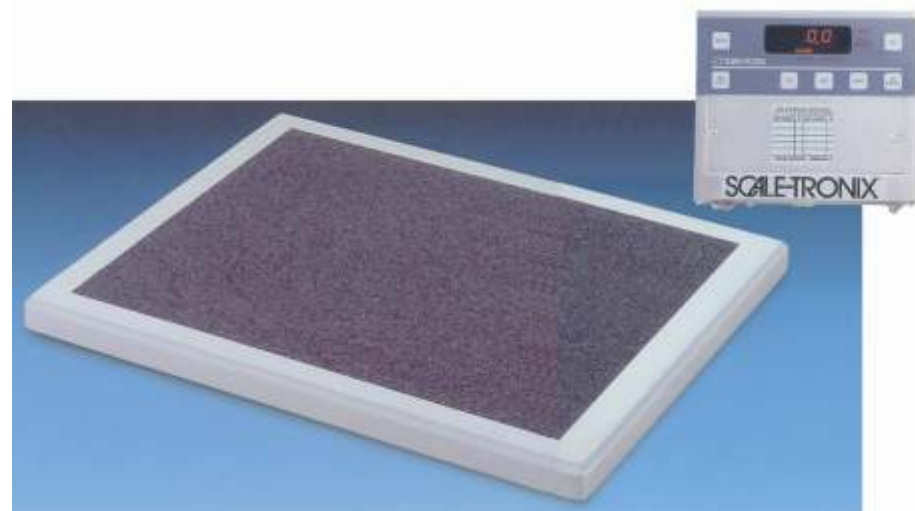

# 5122

## Low Profile Stand-On Scale

After each weighing, the 5122 returns to zero enabling you to weigh a stream of patients without delays. With the press of a button, patients may be reweighed while still on the scale. Accuracy is 1/10 lb (100 gm) with a capacity of 660 lb (300 kg). Power for approximately 20,000 weighings is supplied by 6 "D" alkaline disposable batteries. Line-cord operation is available with an optional power adapter. The new 5122 Stand-On Scale is extremely compact—27" high, 16" wide, 18" deep. It's ideal for use in triage units, clinics, examination and patient rooms. This free-standing digital unit features a 1 1/2" low, non-slip platform and a recessed handle at the rear for easy portability.

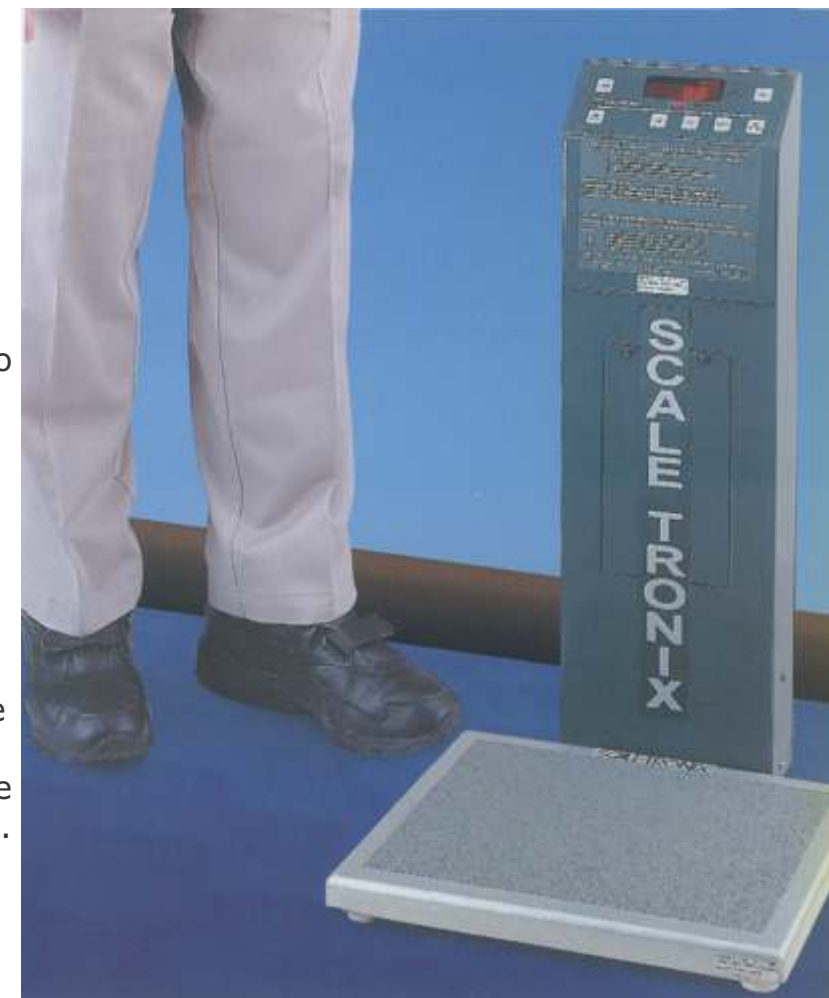

### **Specifications: 5122 Low-Profile Stand-On Scale**

**ACCURACY:** 1/10 LB/100 GM.

**WEIGHING RANGE:** 660 LB/300 KG.

**READOUT:** DIGITAL LED DISPLAY IN POUNDS AND KILOGRAMS.

**POWER SOURCE CORDLESS:** 6 D-SIZE DISPOSABLE ALKALINE BATTERIES. OPTIONAL POWER ADAPTER FOR LINE-CORD OPERATION.

**AUTOMATIC ZERO:** AUTOMATICALLY RETURNS TO ZERO READY FOR NEXT PATIENT

**WEIGHT RECALL:** RECALLS LAST WEIGHT.

**REWEIGH:** RECOMPUTES WEIGHT WHILE PATIENT IS ON THE SCALE

**PLATFORM:** LOW PROFILE 16" WIDE X 14" DEEP X 1 1/2" HIGH

**HEIGHT AND WEIGHT:** 27" HIGH, 19.5 POUNDS.

**LOW BATTERY INDICATOR.**

**WEIGHING MECHANISM:** ELECTRONIC LOAD CELLS

The battery powered 5102 Stand-On Scale is designed for patients up to 660 lbs. It weighs accurately regardless of floor surface, or patient's position on the scale. Weight is displayed on a stainless steel, wall-mounted digital readout. The 1 1/4 inch high platform is stable and is easily moved. A wall bracket provides convenient storage of the platform. Line cord power adapter is optional.

### **Specifications:**

**5102 Physician's Stand-On Scale**

**ACCURACY:** 1/10 LB/100 GM.

**WEIGHING RANGE:** 880 LBS/400 KG.

**READOUT:** DIGITAL LED DISPLAY IN POUNDS AND KILOGRAMS, STAINLESS STEEL, WALL MOUNTED.

**POWER SOURCE:** 6 DISPOSABLE ALKALINE "C" BATTERIES FOR 10,000 WEIGHINGS. EXTERNAL ADAPTER IS AVAILABLE FOR LINE-POWER OPERATION.

**PLATFORM:** 18"w X 14"d X 1 1/4" h.

**WEIGHING MECHANISM:** ELECTRONIC LOAD CELLS.

# 6102 Computer Enhanced Flush-Mounted Floor Scale

- For wheelchair, chair and stand-on weighing.
- Scale is recessed into the floor.
- 825 lb/375 kg capacity.
- Heavy-duty handrail for patient safety.
- Bright, easy-to-read digital display.

The 6102 platform is level with the floor for easy-access, multi-purpose weighing. The patient simply walks, or is wheeled onto the weighing platform and is weighed. A dual system is also available for stretchers, beds, wheelchairs, etc.

Weight is displayed to an accuracy of 100 gm or 1/10 lb on a recessed or surface-mounted digital readout. It's line-cord powered for continual use, with a rechargeable battery backup providing operation during power outages. The scale requires no zeroing between patients.

The optional heavy-duty handrail provides patient safety and convenience, and is positioned to remain clear of wheelchairs and beds.

The 6102 is a computer-enhanced scale providing:

- Digital display in either pounds or kilos.
- Automatic zero.
- Reweigh (recomputes weight while patient is on scale).



## Specifications: 6102 Flush Floor Scale

**ACCURACY:** 1/10 LB, 100 GM. **WEIGHING RANGE:** 825 LB/375 KG.

**READOUT:** DIGITAL DISPLAY IN KILOS AND POUNDS. SURFACE OR FLUSH MOUNTED.

**POWER SOURCE:** 110 VAC, WITH RECHARGEABLE BATTERY BACKUP.

**PLATFORM (DECK):** ANODIZED ALUMINUM, 32" X 36".

6154 Platform (deck): 35" x 54".

**WEIGHING MECHANISM:** ELECTRONIC LOAD CELLS.

**ROUGH-IN FRAME KIT:** ALLOWS FOR A UNIFORM, FINISHED LOOKING INSTALLATION.

**AUTOMATIC ZERO:** RETURNS TO ZERO READY FOR THE NEXT PATIENT.

**REWEIGH:** RECOMPUTES WEIGHT WHILE PATIENT IS ON SCALE.

**WEIGHT RECALL:** RECALLS LAST WEIGHT.

**SUCCESSIVE WEIGHING:** WEIGH A CONTINUOUS STREAM OF PATIENTS WITHOUT RE-ZEROING.

**AUDIBLE SIGNAL:** INDICATES WEIGHT IS IN MEMORY.

**PRINTER:** (ACCESSORY) BUILT-IN, PAPER TAPE, USES STANDARD THERMAL PAPER.

**COMPUTER CAPABILITY:** RS 232 OUTPUT (OPTIONAL).

## Also Available:

### 5152 In Floor Stand In Scale

Same scale with a single 18" x 19½" Platform suitable for ambulatory patients only.

Computer enhanced readout with printer.



Rough In Frame Kit shown below - requires 3 inch depth in floor



## Also available from Global Medical



### Lionville Medication Cart and Point of Care Carts

Lionville's "nurse friendly" electronic locking medication cart is available in three models, 800, 600 or 400, each one designed to meet the requirements & size of the nursing unit.

Customize this cart to your needs by adding patient med drawers, storage drawers, pull out work surfaces, waste bins, sharps containers, cup dispensers..... too many to list.

More than 1800 nursing units across Canada rely on Lionville medication carts to secure and organize their patient's meds.



### 6154 Flush Mounted In Floor Scale for wheelchair, stand on and stretchers

Platform: 35 in. x 54 in.



## Scale-Tronix is the standard for healthcare in Canada.

There are more than 2000 Scale-Tronix scales in use in 400 plus healthcare facilities in Canada. 3 of every 4 customers have purchased multiple scales and more than 90% use more than one model of Scale-Tronix scale.



Scale-Tronix 6002 3 in 1 WheelChair, Stand On & Chair Scale

Also available from Scale-Tronix:  
Scale-Tronix 4802 Pediatric Scale



## GLOBAL MEDICAL PRODUCTS

5230 South Service Rd.  
Burlington, ON L7L 5K2  
info@globalmedical.ca  
www.globalmedical.ca

TEL: (905)634-7799  
TOLL FREE: 1-800-387-6095  
FAX: (905)634-2868  
AU QUÉBEC: 1-800-655-9445

